

# CAREPAD<sup>®</sup>

DRIVING INNOVATION IN SURGERY


PATIENT POSITIONING  
WITH CARE



**CAREPAD<sup>®</sup>**  
DRIVING INNOVATION IN SURGERY

A system

INNOVATIVE COMPLETE VERSATILE

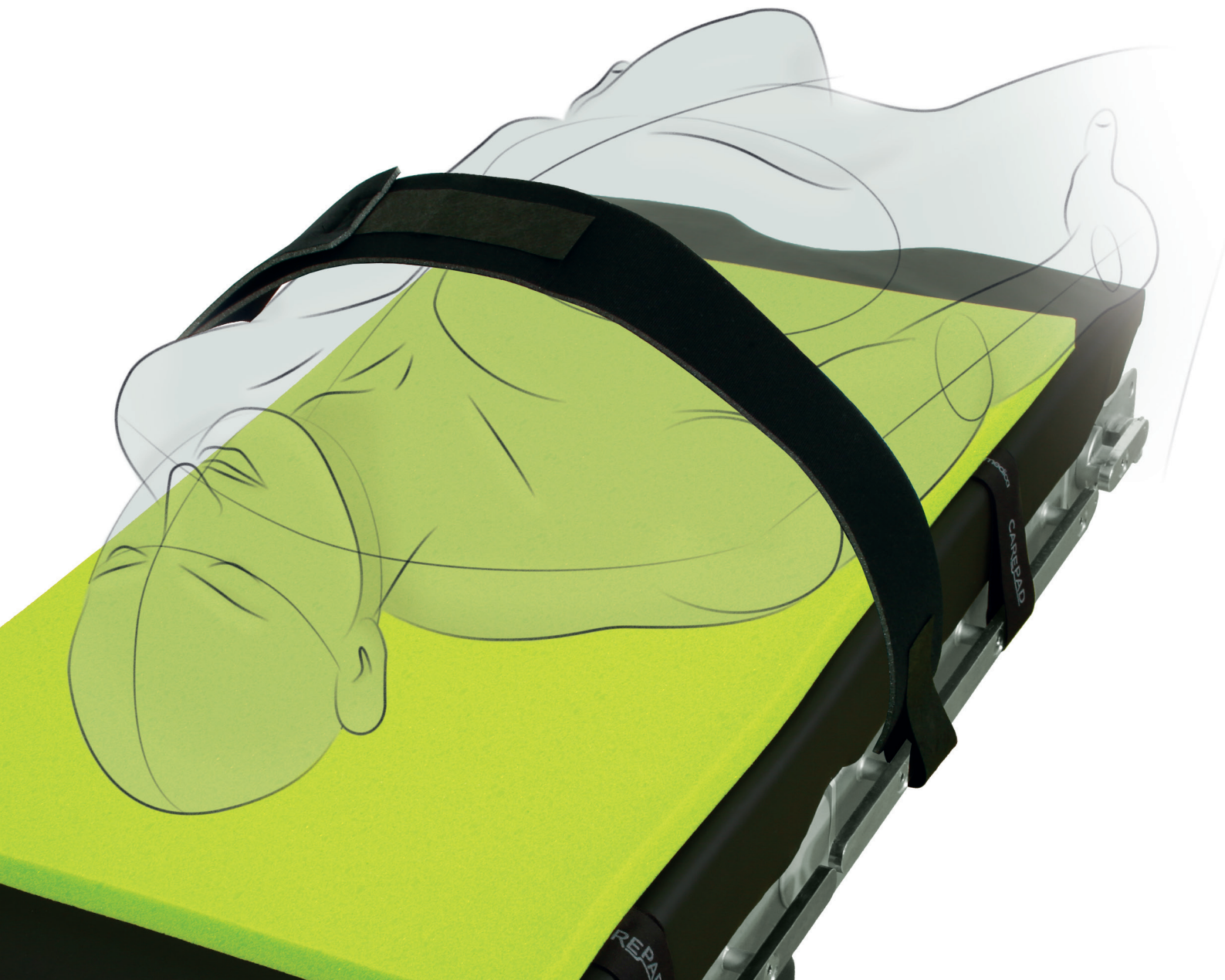



# The system

**CAREPAD®**  
DRIVING INNOVATION IN SURGERY

is an innovative system for positioning the patient during surgical procedures performed with traditional, robotic or laparoscopic techniques.

It allows a Trendelenburg positioning of the patient on the operating bed in a safe, easy and quick manner

## Patient safety

It eliminates the risk of slippage by adapting to the shape and weight of the patient which is uniformly redistributed over the entire mat.

## Efficiency

Simple setup thanks to a patented, repeatable and safe solution.  
The optimal weight redistribution in the pressure zones guarantees excellent performance even without the use of shoulder pads, strips or adhesives.

## Rapidity

Pre-assembled, positionable in a few seconds and ready for use. Quick repositioning thanks to the patented solution.

## Biocompatibility<sup>1</sup>


Biocompatible<sup>1</sup> and free of latex or phthalates<sup>2</sup>, it is made of viscoelastic Memory material with slow return and high density, offering thermal insulation and anti-decubitus function.

## Trendelenburg positioning

CarePad is ideal for the patient's in Trendelenburg positioning pushed up to 45 degrees of inclination.

## Medical Device

CarePad is a Medical Device, CE marked, in accordance with Regulation (EU) 2017/745

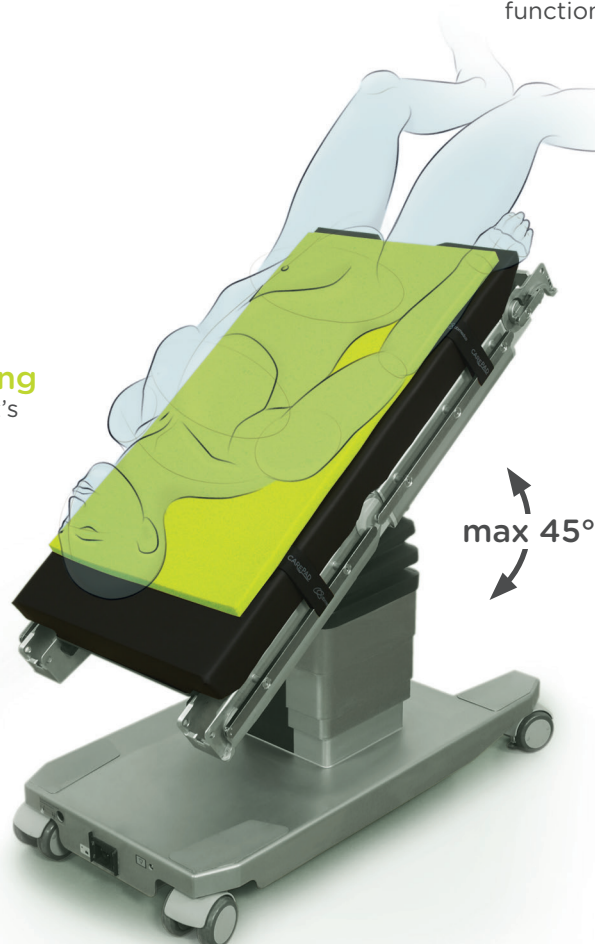

<sup>1</sup> Biocompatibility assessed on the device, in accordance with ISO 10993-1:2018

<sup>2</sup> No substances containing latex and its derivatives or phthalates are intentionally used in its manufacture



# The material

The CarePad mat, biocompatible<sup>1</sup> and free of latex or phthalates<sup>2</sup>, is made of slow-return, high-density viscoelastic Memory material.

Its particular composition with open cells guarantees thermal insulation while maintaining the patient's body temperature constant.

Its anti-decubitus function reduces the risk of post-operative redness and irritation and pressure ulcers.


## Viscoelastic material

that adapts to the shape and weight of the patient in order to prevent slippage.



## Patented repositioning

Thanks to a patented solution, a simple repositioning of the velcro straps over the entire length of the mattress is allowed to enable a perfect and customized locking to the bed bars.

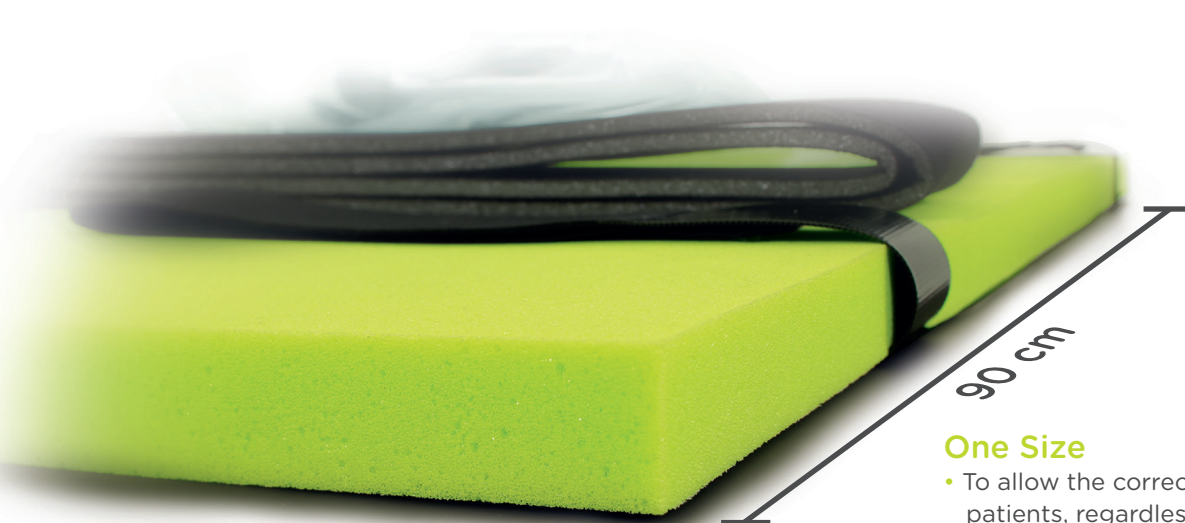



## One Size

- To allow the correct positioning of patients, regardless of height and weight (maximum weight 250 Kg).
- To facilitate orders and storage thanks to the use of a single code.

# The kit

Designed to be as complete and versatile as possible, it includes:

- CarePad mat with repositionable velcro side straps
- 2 additional velcro straps
- 1 patient sheet
- 1 chest strap

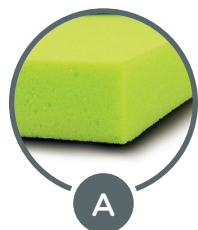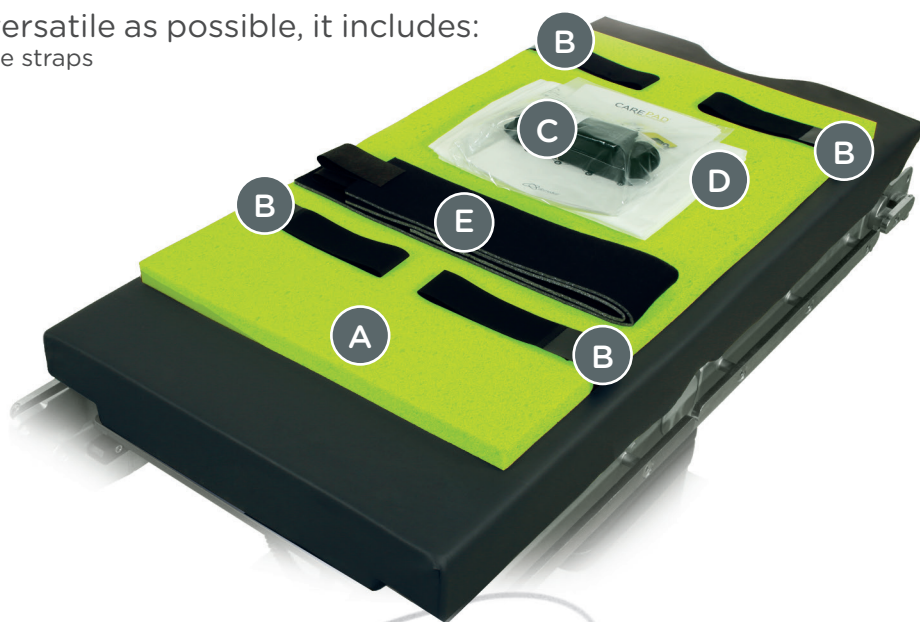

## Mat

Thanks to its length of 90 cm, it allows the correct positioning of patients regardless of height and weight (max weight 250 Kg).



## Side straps

Easy and versatile anchoring system to the operating bed.

High-quality velcro, pre-positioned on the mat, reconfigurable in case of need by the operator thanks to the patented solution.



## Additional belts

Additional longer length straps to make the mat easily positionable on operator's beds without side bars or for extra fixing.



## Patient sheet

Sheet for securing the patient's arms and for repositioning the patient on the mat.




## Thoracic strap

Provides additional patient safety.

With no buttonholes, it is extremely easy and quick to install.





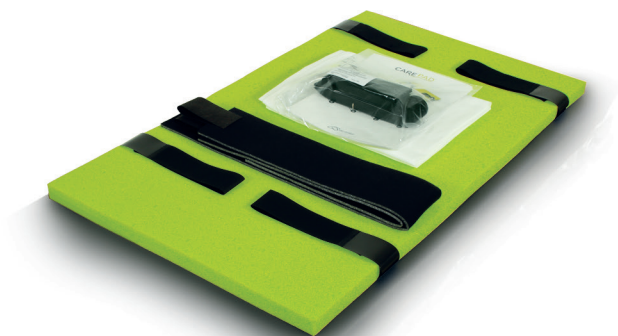

**CAREPAD<sup>®</sup>**  
DRIVING INNOVATION IN SURGERY

PATENT N°

IT 102020000006058 / IT 102020000024568 / IT 202020000001249 / IT 202020000005827

#### Codes number

Code	Description	Piece per pack
3/PAD90503	CAREPAD Kit for patient's Trendelenburg positioning	6
3/PAD102503	CAREPAD XL Kit for patient's Trendelenburg positioning	3

#### The kit includes:

- 1 **CAREPAD** mat with repositionable velcro side straps  
Dimensions: 90x50x2.5 cm
- 2 additional velcro straps
- 1 patient sheet
- 1 chest strap



**ab medica s.p.a.**

Via J.F. Kennedy, 10/12 - Cerro Maggiore (MI) - Italy  
tel. +39 02933051 - fax +39 0293305400  
[www.abmedica.it](http://www.abmedica.it)