

# LiNA MaxFlow™ & LiNA MaxFlow™ Endohook

Powerful suction and irrigation system  
with electrosurgery applications



# LiNA MaxFlow™ & LiNA MaxFlow™ Endohook

**Safe.** Designed for suction and irrigation.

Suction probe with relieving eyes and false-air regulator on suction handle designed for procedures where delicate suction is required

LiNA MaxFlow™ Endohook

Monopolar endohook with safety function for retraction of electrode tip into cannula.

**Versatile.** Variety of lengths, electrode tips and surgical applications.

Suction probe available in 5mm or 10mm, and 5mm extended length (370mm)

Choice of solely suction and irrigation function or supplement with electrosurgical monopolar accessory provided by LiNA MaxFlow™ Endohook

LiNA MaxFlow™ Endohook

3 different electrode tips on monopolar endohook



Safety on/ready  
for suction  
and irrigation



Safety off/ready  
for electrosurgery


Click-on locking mechanism




LiNA MaxFlow™  
Endohook



Relieving eyes

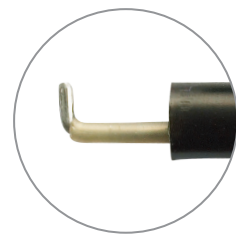


L-shaped tip

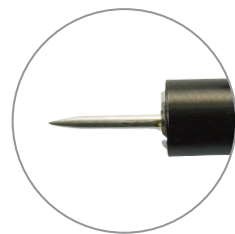

Needle tip



Spatula tip

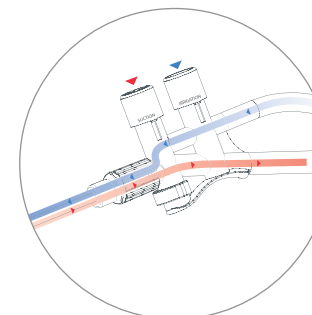

**Functional.** Beneficial features to meet surgical needs.

Unique separate channel design in handle designed to reduce risk of blood clots or debris blockage


True 9mm direct flow through whole suction system when using a 10mm probe during suction and irrigation procedure

No need for frequent change of instruments during surgical procedure as suction, irrigation, dissection, coagulation and cutting is all in one instrument!



# LiNA MaxFlow™

## & LiNA MaxFlow™ Endohook

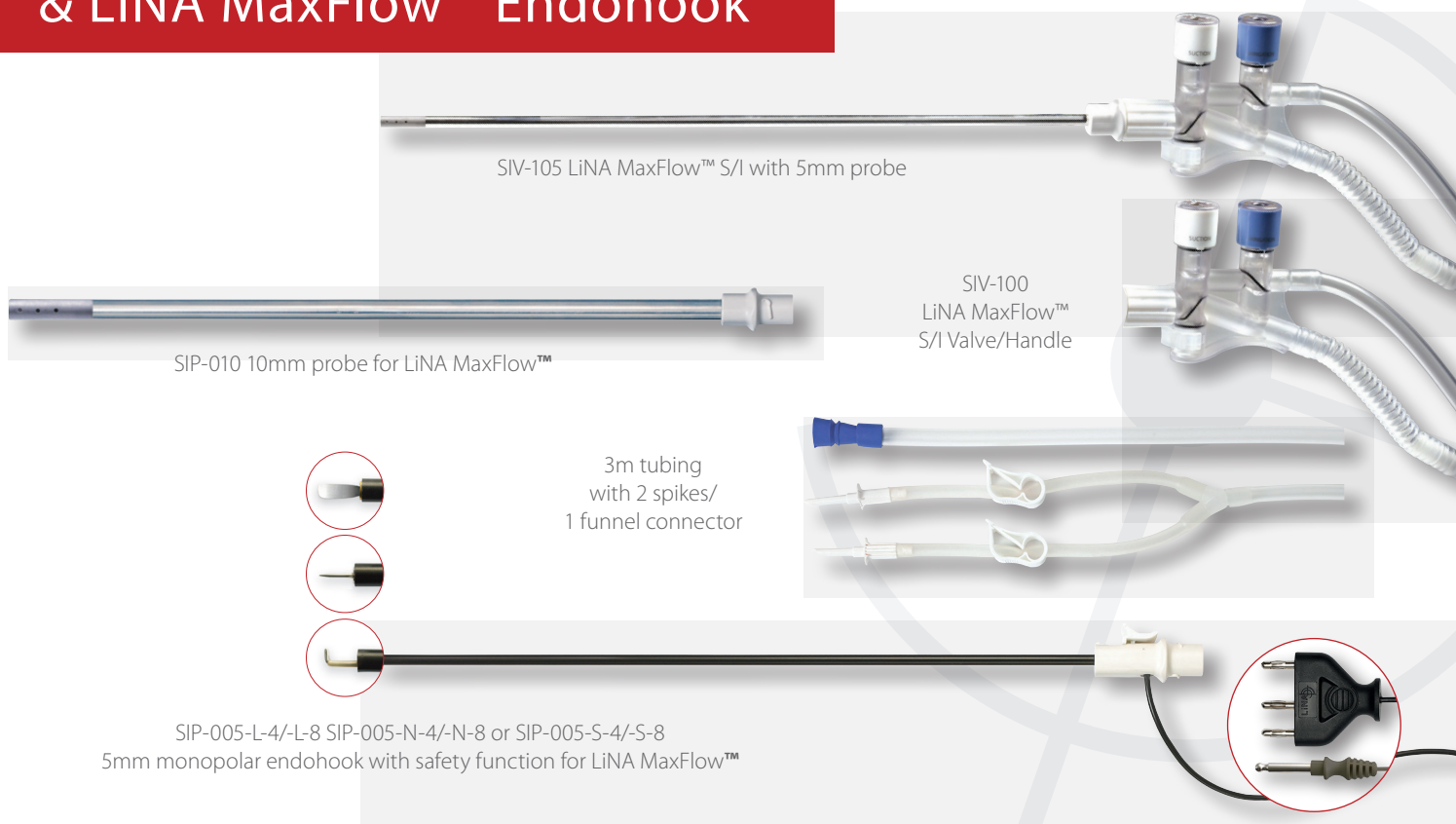

### Ordering information:

LiNA MaxFlow™ **CE**  
0459

| Product code | Description   | UOM          |
|--------------|---|--------------|
| SIV-100      | LiNA MaxFlow™ (valve only) – 3m tubing and 2 spikes plus connector. Sterile, single use.  | Box of 10 ea |
| SIV-105      | LiNA MaxFlow™ suction/irrigation set – with 5mm probe, 3m tubing and 2 spikes plus connector. Sterile, single use.                            | Box of 10 ea |
| SIP-005-370  | 5 mm probe, extended length (370 mm). Sterile, single use.  | Box of 10 ea |
| SIP-010      | 10 mm probe. Sterile, single use.   | Box of 10 ea |
| SIP-005-L-4  | LiNA MaxFlow™ Endohook – Endoscopic monopolar hook for LiNA MaxFlow™, 3m cable, 3-pin connector, L-shaped tip electrode. Sterile, single use. | Box of 10 ea |
| SIP-005-L-8  | LiNA MaxFlow™ Endohook – Endoscopic monopolar hook for LiNA MaxFlow™, 3m cable, 8 mm connector, L-shaped tip electrode. Sterile, single use.  | Box of 10 ea |
| SIP-005-S-4  | LiNA MaxFlow™ Endohook – Endoscopic monopolar hook for LiNA MaxFlow™, 3m cable, 3-pin connector, Spatula tip electrode. Sterile, single use.  | Box of 10 ea |
| SIP-005-S-8  | LiNA MaxFlow™ Endohook – Endoscopic monopolar hook for LiNA MaxFlow™, 3m cable, 8 mm connector, Spatula tip electrode. Sterile, single use.   | Box of 10 ea |
| SIP-005-N-4  | LiNA MaxFlow™ Endohook – Endoscopic monopolar hook for LiNA MaxFlow™, 3m cable, 3-pin connector, Needle tip electrode. Sterile, single use.   | Box of 10 ea |
| SIP-005-N-8  | LiNA MaxFlow™ Endohook – Endoscopic monopolar hook for LiNA MaxFlow™, 3m cable, 8 mm connector, Needle tip electrode. Sterile, single use.    | Box of 10 ea |

### Distributor:

©2022 LiNA Medical ApS. All rights reserved.

LiNA is a registered trademark of LiNA Medical in EU, US and CN.

LiNA Medical ApS

• Formervangen 5  
DK-2600 Glostrup  
Denmark

• Tel.: +45 43 29 66 66  
Fax: +45 43 29 66 99  
info@lina-medical.com

• www.lina-medical.com  
CVR: 10054974

EN ISO 13485:2016  
Quality Management  
System

