

LiNA Port™

Single use 5mm trocar system



LiNA Port™

Single use 5mm trocar system

The LiNA Port™ is a single use trocar system fulfilling the demand for quality and efficiency at a competitive price.



Safety. **Silicone valve – airtight seal**
Ensuring minimum loss of pneumoperitoneum during surgery.

Dilating tip for delicate tissue
The blunt dilating tip in the LiNA Port™ range supports non-traumatic entry.

Functionality. **Threaded cannula**
To optimize efficient abdominal retention.



Versatility. **Choice of Dilating, Conical or Pyramidal Tip**
Beside the dilating tip, the LiNA Port™ is available with conical or pyramidal tip for sharp and easy entry.

Two Lengths - Short and Long Ports
Short 60 mm and Long 87 mm. Internal cannula diameter: ø5,8 mm allows easy insertion of all 5mm instruments.
Suitable for laparoscopic procedures on most patients.

Gas Ports
The cannulas are available with either a simple gas port, a gas port with on/off valve or without gas port.


Single Pack and Dual Pack
The LiNA Port™ is available in a single or dual pack. The dual pack includes one trocar, two cannulas and one LiNA Clear Sight Wipe (LiNA AntiFog Solution™) for camera optics.

Airtight Silicone Double Seal Valve.
Designed for secure position of instrument, easy exchange of instruments and enables the use of instruments down to 3mm.



Non-traumatic Fascia Thread.
Designed for safe and firm anchoring of cannula.

Dilating, Conical or Pyramidal Tip.
Supports the most used entry techniques.

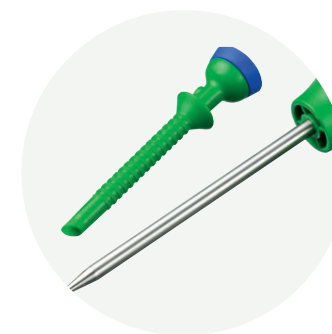

Single Pack
LiNA Port™
5mm trocar system



Dual Pack
with one
LiNA Clear Sight Wipe
(LiNA AntiFog Solution™)



LiNA Port™

Single use 5mm trocar system

LiNA Port™  0459

Ordering Information:

REF.	Dilating Tip	UOM
1005-06D1	5mm, short, dilating trocar + w/o gas port cannula. Sterile, single use.	Box of 5 sets
1005-06D2	5mm, short, dilating trocar + 2 x cannula w/o gas port + wipe. Sterile, single use.	Box of 5 sets
1005-09D1	5mm, long, dilating trocar + cannula w/o gas port. Sterile, single use.	Box of 5 sets
1005-09D2	5mm, long, dilating trocar + 2 x cannula w/o gas port + wipe. Sterile, single use.	Box of 5 sets
1005-06D1-GV	5mm, short, dilating trocar + cannula w/valve. Sterile, single use.	Box of 5 sets
1005-06D2-GV	5mm, short, dilating trocar + cannula w/valve + cannula w/o gas port + wipe. Sterile, single use.	Box of 5 sets
1005-09D1-GV	5mm, long, dilating trocar + cannula w/valve. Sterile, single use.	Box of 5 sets
1005-09D2-GV	5mm, long, dilating trocar + cannula w/valve + cannula w/o gas port + wipe. Sterile, single use.	Box of 5 sets
1005-06D1-G	5mm, short, dilating trocar + cannula w/gas port. Sterile, single use.	Box of 5 sets
1005-06D2-G	5mm, short, dilating trocar + cannula w/gas port + cannula w/o gas port + wipe. Sterile, single use.	Box of 5 sets
1005-09D1-G	5mm, long, dilating trocar + cannula w/gas port. Sterile, single use.	Box of 5 sets
1005-09D2-G	5mm, long, dilating trocar + cannula w/gas port + cannula w/o gas port + wipe. Sterile, single use.	Box of 5 sets

REF.	Conical Tip	UOM
1007-1	5mm, short, conical trocar + cannula w/o gas port. Sterile, single use.	Box of 5 sets
1007-2	5mm, short, conical trocar + 2 x cannula w/o gas port + wipe. Sterile, single use.	Box of 5 sets
1007-2-2GP	5mm, short, conical trocar + 2 x cannula w/valve + wipe. Sterile, single use.	Box of 5 sets
1009-C-1	5mm, short, conical trocar + cannula w/gas port. Sterile, single use.	Box of 5 sets
1009-C-2GP	5mm, short, conical trocar + 2 x cannula w/gas port + wipe. Sterile, single use.	Box of 5 sets
2007-1	5mm, long, conical trocar + cannula w/o gas port. Sterile, single use.	Box of 5 sets
2007-2	5mm, long, conical trocar + 2 x cannula w/o gas port + wipe. Sterile, single use.	Box of 5 sets

REF.	Pyramidal Tip	UOM
1005-1	5mm, short, pyramidal trocar + cannula w/o gas port. Sterile, single use.	Box of 5 sets
1005-2	5mm, short, pyramidal trocar + 2 x cannula w/o gas port + wipe. Sterile, single use.	Box of 5 sets
2005-1	5mm, long, pyramidal trocar + cannula w/o gas port. Sterile, single use.	Box of 5 sets
2005-2	5mm, long, pyramidal trocar + 2 x cannula w/o gas port + wipe. Sterile, single use.	Box of 5 sets
1008-1	5mm, short, pyramidal trocar + cannula w/valve. Sterile, single use.	Box of 5 sets
1008-2	5mm, short, pyramidal trocar + cannula w/valve + cannula w/o gas port + wipe. Sterile, single use.	Box of 5 sets
2008-1	5mm, long, pyramidal trocar + cannula w/valve. Sterile, single use.	Box of 5 sets
2008-2	5mm, long, pyramidal trocar + cannula w/valve + cannula w/o gas port + wipe. Sterile, single use.	Box of 5 sets
1009-1	5mm, short, pyramidal trocar + cannula w/gas port. Sterile, single use.	Box of 5 sets
1009-2	5mm, short, pyramidal trocar + cannula w/gas port + cannula w/o gas port + wipe. Sterile, single use.	Box of 5 sets
2009-1	5mm, long, pyramidal trocar + cannula w/gas port. Sterile, single use.	Box of 5 sets
2009-2	5mm, long, pyramidal trocar + cannula w/gas port + cannula w/o gas port + wipe. Sterile, single use.	Box of 5 sets

Distributor:

©2020 LiNA Medical ApS. All rights reserved.
LiNA is a registered trademark of LiNA Medical in EU, US and CN.

LiNA Medical ApS

• Formervangen 5
DK-2600 Glostrup
Denmark

• Tel.: +45 43 29 66 66
Fax: +45 43 29 66 99
info@lina-medical.com

• www.lina-medical.com
CVR: 10054974

Quality Management System
EN ISO13485:2016

